

ASAHI

PTCA Guide Wire

ASAHI's unique guide wire technology with Composite Core

SION TECC 

Leading to the NEXT

ASAHI

PTCA Guide Wire

The Right Wire for Any Case


The reliability of the products is highly recognized by its users in the medical profession.

SION TECC

ASAHI's unique guide wire technology with Composite Core




- **Enhanced tip durability**
- **Combination of flexibility and torque force**
- **Reduced whipping of the wire due to superior torque response**

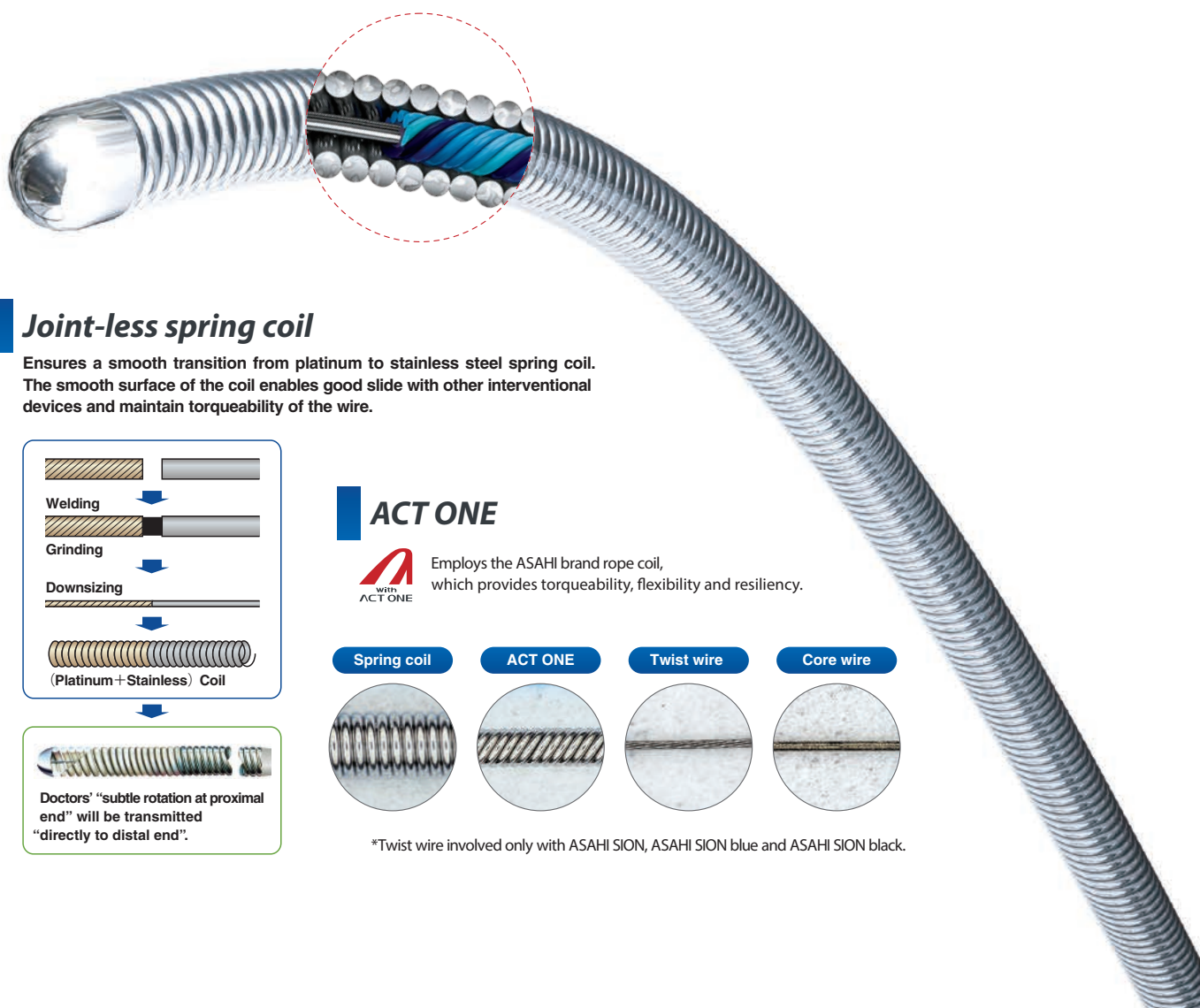

Joint-less spring coil

Ensures a smooth transition from platinum to stainless steel spring coil. The smooth surface of the coil enables good slide with other interventional devices and maintain torqueability of the wire.

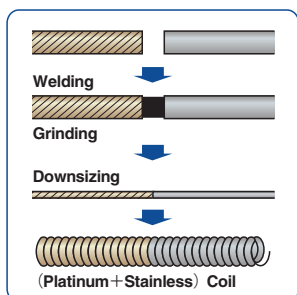


Doctors' "subtle rotation at proximal end" will be transmitted "directly to distal end".

ACT ONE



Employs the ASAHI brand rope coil, which provides torqueability, flexibility and resiliency.

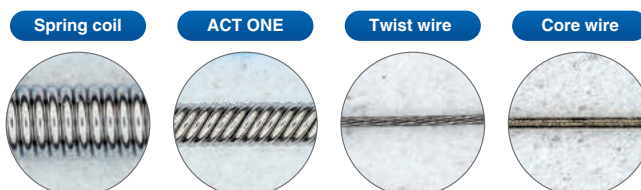

*Twist wire involved only with ASAHI SION, ASAHI SION blue and ASAHI SION black.

Tip design


A variety of tip designs to facilitate branch selection or penetration ability. Asahi offers the right wire for any lesion.

One-piece core wire

One-piece core wire supports the entire guide wire from the proximal to the distal end. This design transmits the guide wire torque from one end to the other without loss.


Coatings

A variety of coatings provide either superior lubricity, tactile feedback, or an optimal combination of the two.

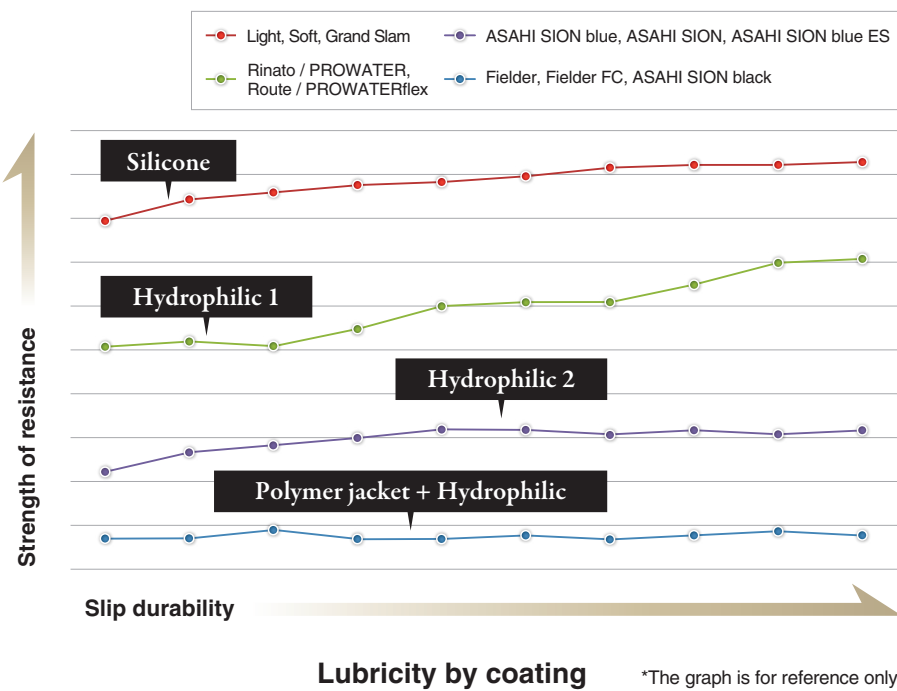

Complex and Chronic occlusion guide wires with SION TECC


ASAHI SUOH 03

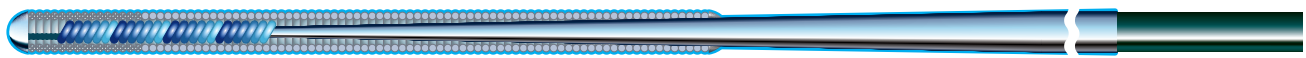

- Tip load 0.3 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating 52 cm

Flexible shaft provides superb trackability and crossability in severe tortuosity. The long radiopaque segment provides good visualization.

Fielder XT-R

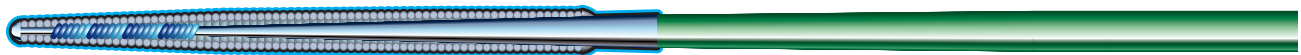

- Tip load 0.6 gf
- Tip radiopacity 16 cm
- Polymer jacket length 17 cm
- Tip outer diameter 0.26 mm (0.010 inch)
- SLIP-COAT® coating over the polymer jacket

Guide wire recommended for narrow channel tracking with its low profile, flexible tip, and high lubricity polymer coating.

Fielder XT-A



- Tip load 1.0 gf
- Tip radiopacity 16 cm
- Polymer jacket length 17 cm
- Tip outer diameter 0.26 mm (0.010 inch)
- SLIP-COAT® coating over the polymer jacket

Guide wire with a higher tip load than Fielder XT and XT-R, facilitating entry into the chronic occluded lesion.

ASAHI Gaia series (First, Second, Third)




- Tip load
 - ASAHI Gaia First 1.7 gf
 - ASAHI Gaia Second 3.5 gf
 - ASAHI Gaia Third 4.5 gf
- Tip radiopacity 15 cm
- Tip outer diameter
 - ASAHI Gaia First 0.26 mm (0.010 inch)
 - ASAHI Gaia Second 0.28 mm (0.011 inch)
 - ASAHI Gaia Third 0.30 mm (0.012 inch)
- SLIP-COAT® coating 40 cm

Guide wire with a high maneuverability in the chronic total occluded lesion, and a tip designed to improve penetrability into the lesion while remaining flexible.

Frontline guide wires


Rinato / PROWATER

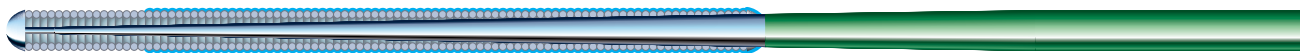

- Tip load 0.8 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating over the spring coil*

First choice guide wire with floppy tip and moderate core wire support. Excellent maneuverability and good support for device delivery.
*Silicone coating up to 30mm from the tip for safer procedures.


Route / PROWATERflex

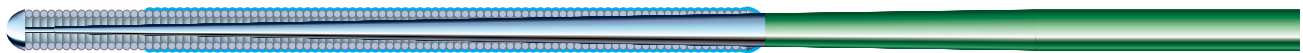

- Tip load 0.8 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating over the spring coil*

First choice guide wire featuring a soft tip with good shaping memory and excellent guide wire trackability. Sufficient support for delivery of most interventional devices.
*Silicone coating up to 30mm from the tip for safer procedures.


Soft

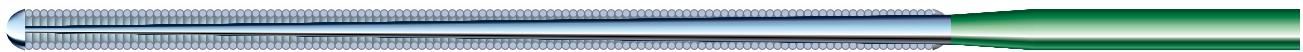

- Tip load 0.7 gf
- Tip radiopacity 3 cm
- Silicone coating over the spring coil

Basic frontline guide wire with floppy tip and non hydrophilic coating to ensure optimal tactile feedback.

Light

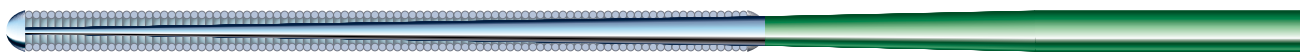

- Tip load 0.5 gf
- Tip radiopacity 3 cm
- Silicone coating over the spring coil

Basic frontline guide wire with a flexible core and good torqueability for distal side branch selection.

Felder



- Tip load 1.0 gf
- Tip radiopacity 3 cm
- Polymer jacket length 22 cm
- SLIP-COAT® coating over the polymer jacket

A polymer jacket guide wire, providing advanced slip performance and trackability.

Felder FC



- Tip load 0.8 gf
- Tip radiopacity 3 cm
- Polymer jacket length 20 cm
- SLIP-COAT® coating over the polymer jacket

Fine control over challenging tortuous vessels and highly stenosed lesions. Polymer jacket provides advanced slip performance with superior torque and support.

Others

Grand Slam



- Tip load 0.7 gf
- Tip radiopacity 4 cm
- Silicone coating over the spring coil

Stiff wire shaft provides extra support for delivery of interventional devices in extreme cases.

Complex and Chronic occlusion guide wires

Fielder XT



- Tip load 0.8 gf
- Tip radiopacity 16 cm
- Polymer jacket length 16 cm
- Tip outer diameter 0.23 mm (0.009 inch)
- SLIP-COAT® coating over the polymer jacket

Polymer jacket, providing excellent lubricity and trackability in tortuous vessels. The tapered tip provides extreme precision for the treatment of complex lesions such as sub-total occlusions and long diffused lesions.

Intermediate / MEDIUM



- Tip load 3.0 gf
- Tip radiopacity 3 cm
- Silicone coating over the spring coil

First choice wire for chronic occlusion cases.

ULTIMATEbros 3



- Tip load 3.0 gf
- Tip radiopacity 11 cm
- SLIP-COAT® coating over the spring coil until the middle of the shaft / without tip

Long hydrophilic coating maintains high maneuverability, allowing for improved wire manipulation in heavy stenosed lesions. Fine shaping improves vessel selectivity and reduces the risk of false lumen expansion.

ASAHI Miracle series / ASAHI MIRACLEbros series

3 / 4.5 / 6 / 12



- Tip load 3.0 / 4.5 / 6.0 / 12.0 gf
- Tip radiopacity 11 cm
- Silicone coating over the spring coil

Guide wire with good tactile feed back which maintains torqueability through the occlusion.

ASAHI Conquest / ASAHI Confianza



- Tip load 9.0 gf
- Tip radiopacity 20 cm
- Silicone coating over the spring coil

Silicone coated and tapered to 0.23mm (0.009 inch) at the tip. The 9gf stiff, tapered tip helps to penetrate highly stenosed lesions.

ASAHI Conquest Pro series / ASAHI Confianza Pro series

9 / 12 / 8-20



- Tip load 9.0 / 12.0 / 20.0 gf
- Tip radiopacity 20 cm (8-20 17 cm)
- SLIP-COAT® coating over the spring coil, excluding the tip
- Tip outer diameter
ASAHI Conquest Pro, Pro12 0.23 mm (0.009 inch)
ASAHI Conquest Pro 8-20 0.20 mm (0.008 inch)

Tapered guide wire for crossing complex lesion with heavy calcification or tough fibrous tissues.

Externalization guide wire

ASAHI RG3




- Tip load 3.0 gf
- Tip radiopacity 3 cm
- SLIP-COAT® coating over the spring coil and until the middle of the shaft
- Shaft diameter 0.26 mm (0.010 inch)
- Guide wire length 330 cm

For guide wire externalization only. Optimal wire strength, hydrophilic coating and 0.26 mm (0.010 inch) shaft provide superior inside-catheter pushability. With the inner wall damage possibility reduced in tortuous vessels as well, the risk of complication is minimized.

ASAHI's unique guide wire technology with Composite Core


SION TECC

Leading to the NEXT

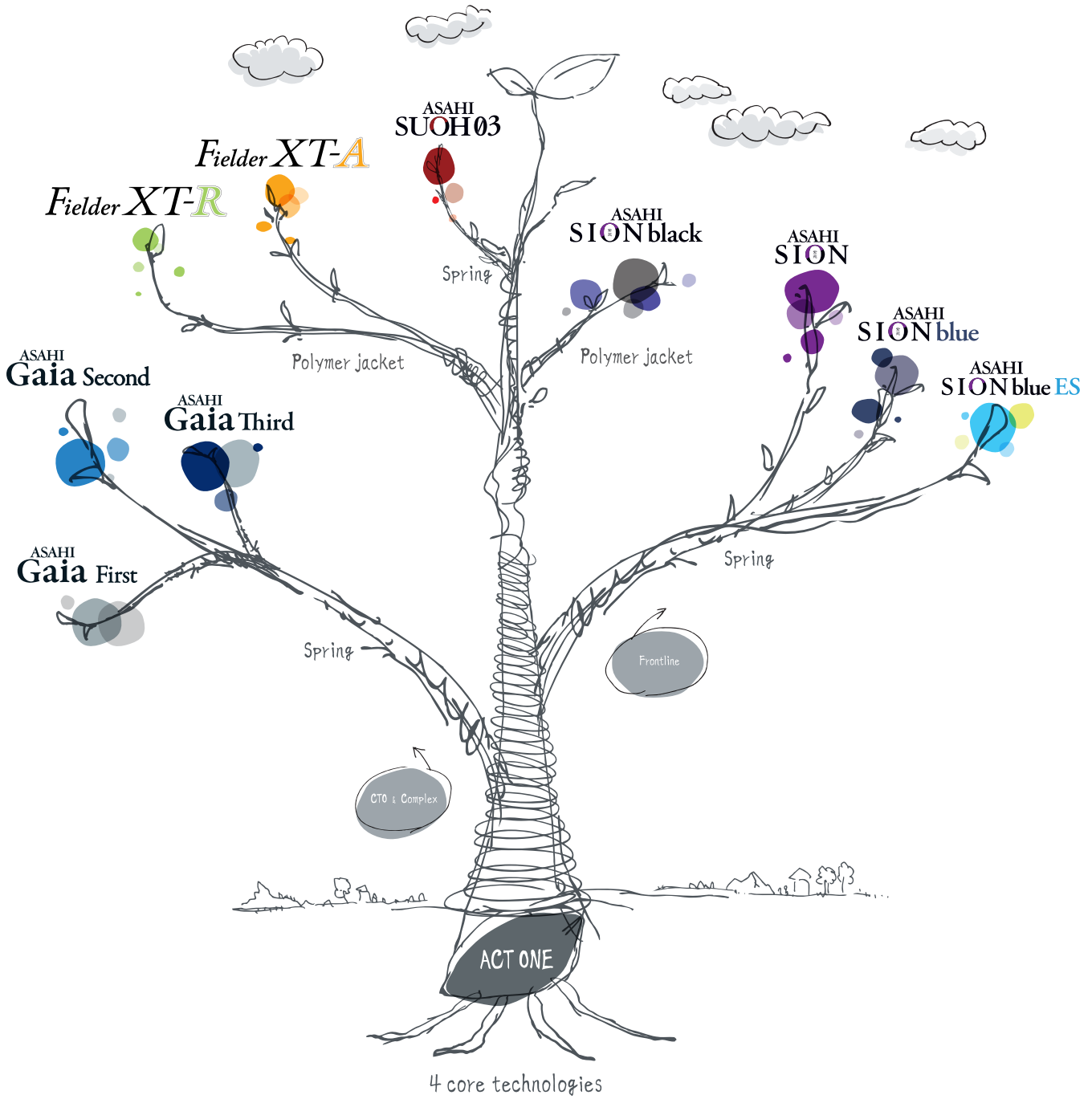


 ASAHI INTECC


Frontline wires

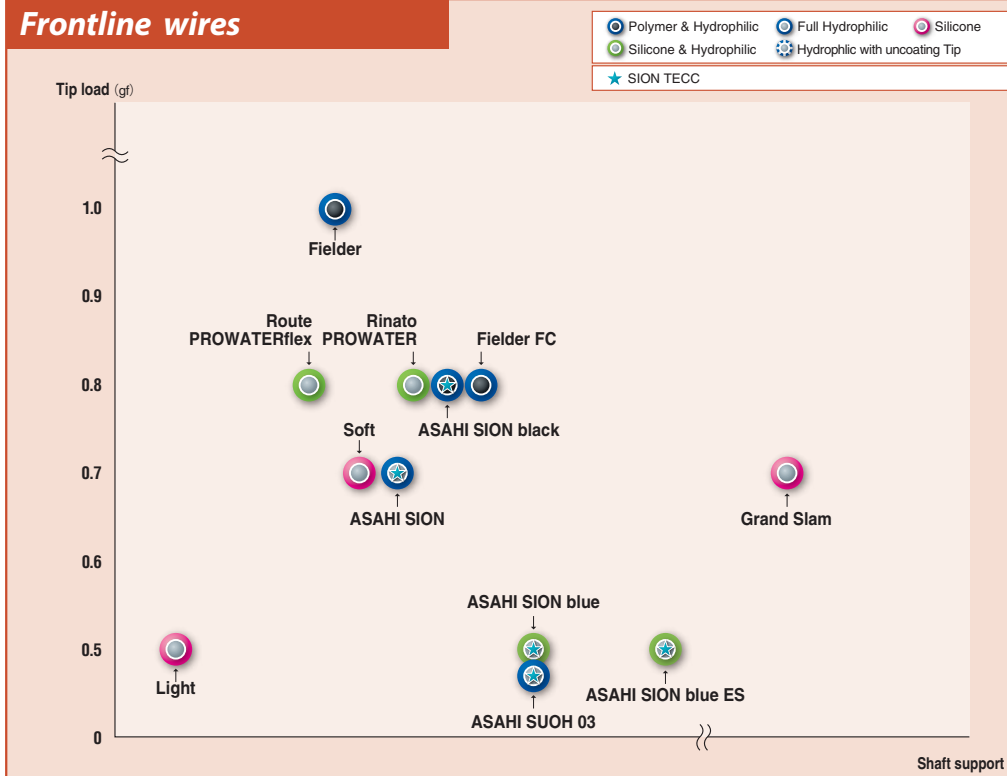

*ASAHI SUOH 03 is recommended to use mainly for the retrograde wiring.

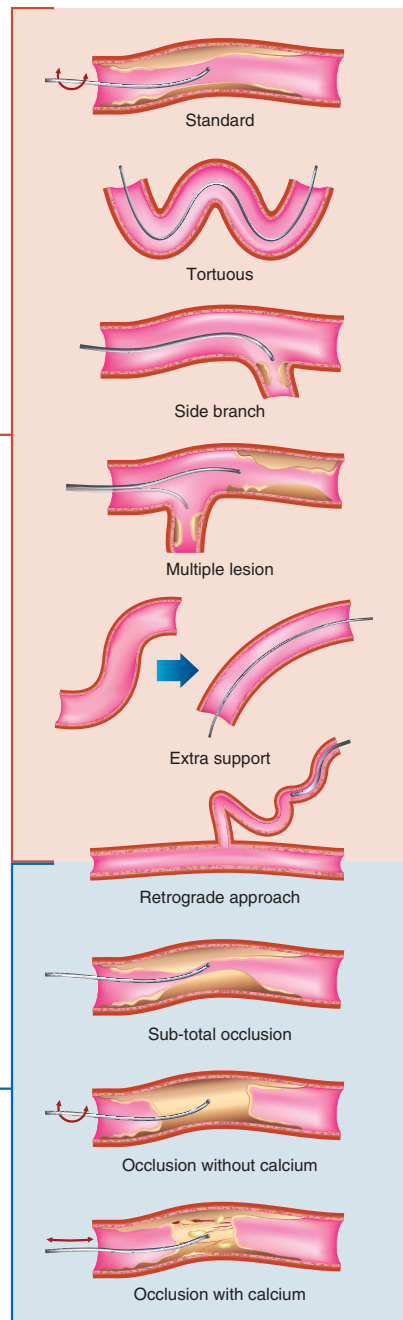

Chronic occlusion wires

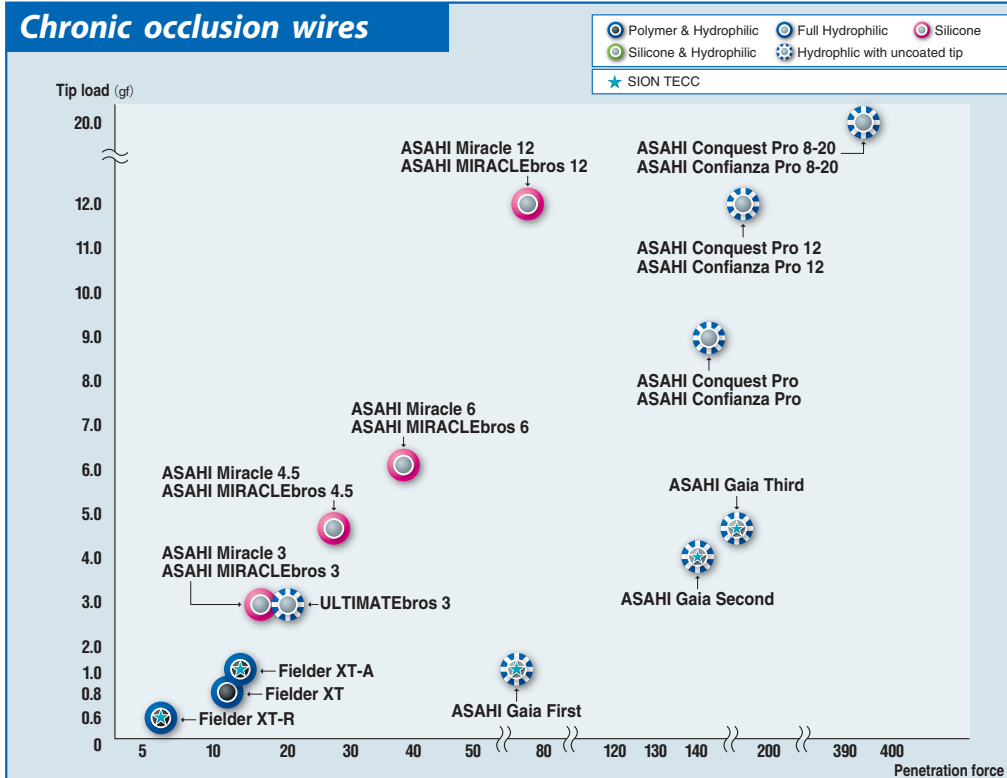

Polymer & Hydrophilic
Enhanced lubricity for better crossability of severe stenoses.

Silicone
Prevent slippage to gain tactile response for lesion morphology.

Silicone & Hydrophilic
Combine advantages of both. Silicone coating over the last few centimeters distally, hydrophilic coating over the rest of the coil and shaft.

Full Hydrophilic
Hydrophilic with uncoated tip
Enhance lubricity for easier wire manipulation in lesions.

SION TECC

Ordering information

SION TECC series

Product	Catalog No.	Diameter (mm (inch))	Radiopaque segment (cm)	Spring coil length (cm)	Length (cm)	Tip shape
ASAHI SION blue	AHW14R004S	0.36 (0.014)	3	20	180	Straight
	AHW14R304S				300	Straight
	AHW14R004J				180	J
	AHW14R304J				300	J
ASAHI SION	AHW14R001S	0.36 (0.014)	3	28	180	Straight
	AHW14R301S				300	Straight
	AHW14R001J				180	J
	AHW14R301J				300	J
ASAHI SION black	APW14R010S	0.36 (0.014)	3	12	180	Straight
	APW14R310S				300	Straight
	APW14R010J				180	J
	APW14R310J				300	J
ASAHI SION blue ES	AHW14R017S	0.36 (0.014)	3	8.5	180	Straight
	AHW14R317S				300	Straight
	AHW14R017J				180	J
	AHW14R317J				300	J
Fielder XT-A	APW14R009S	0.26 / 0.36	16	16	190	Straight
	APW14R309S	(0.010 / 0.014)			300	Straight
Fielder XT-R	APW14R005S	0.26 / 0.36	16	16	190	Straight
	APW14R305S	(0.010 / 0.014)			300	Straight
ASAHI Gaia First	AHW14R007P	0.26 / 0.36	15	15	190	Pre-shape
	AHW14R307P	(0.010 / 0.014)			300	Pre-shape
ASAHI Gaia Second	AHW14R008P	0.28 / 0.36	15	15	190	Pre-shape
	AHW14R308P	(0.011 / 0.014)			300	Pre-shape
ASAHI Gaia Third	AHW14R011P	0.30 / 0.36	15	15	190	Pre-shape
	AHW14R311P	(0.012 / 0.014)			300	Pre-shape
ASAHI SUOH 03	AHW14R013P	0.36 (0.014)	3	19	190	Pre-shape
	AHW14R313P				300	Pre-shape
	AHW14R013S				190	Straight
	AHW14R313S				300	Straight

Traditional wires

Product	Catalog No.	Diameter (mm (inch))	Radiopaque segment (cm)	Spring coil length (cm)	Length (cm)	Tip shape
Route / PROWATERflex	AGH147000	0.36 (0.014)	3	20	180	Straight
	AGH146000	0.36 (0.014)	3	20	180	Straight
Rinato / PROWATER	AGH146300				300	Straight
	AGH146000J	180	J			
Soft	AGH146000	0.36 (0.014)	3	30	180	Straight
	AGH146300				300	Straight
	AGH146000J				180	J
Light	AG145000	0.36 (0.014)	3	20	180	Straight
	AG145300				300	Straight
	AG145000J				180	J
Fielder	AGP140000	0.36 (0.014)	3	12	180	Straight
Fielder FC	AGP140001	0.36 (0.014)	3	11	180	Straight
	AGP140301				300	Straight
Fielder XT	AGP140002	0.23 / 0.36	16	16	190	Straight
	AGP140302	(0.009 / 0.014)			300	Straight
	AG142000	180			Straight	
Intermediate / MEDIUM	AG142000	0.36 (0.014)	3	30	180	Straight
	AG142000J				300	Straight
	AG142000J				180	J
ASAHI Miracle 3 / ASAHI MIRACLEbros 3	AG14M050	0.36 (0.014)	11	11	180	Straight
ASAHI Miracle 4.5 / ASAHI MIRACLEbros 4.5	AG14M045	0.36 (0.014)	11	11	180	Straight
ASAHI Miracle 6 / ASAHI MIRACLEbros 6	AG14M060	0.36 (0.014)	11	11	180	Straight
ASAHI Miracle 12 / ASAHI MIRACLEbros 12	AG14M070	0.36 (0.014)	11	11	180	Straight
ULTIMATEbros 3	AHW14S003S	0.36 (0.014)	11	11	180	Straight
ASAHI Conquest / ASAHI Confianza	AG143090	0.23 / 0.36 (0.009 / 0.014)	20	20	180	Straight
ASAHI Conquest Pro / ASAHI Confianza Pro	AGH143090	0.23 / 0.36 (0.009 / 0.014)	20	20	180	Straight
ASAHI Conquest Pro 12 / ASAHI Confianza Pro 12	AGH143091	0.23 / 0.36 (0.009 / 0.014)	20	20	180	Straight
*ASAHI Conquest Pro 8-20 / ASAHI Confianza Pro 8-20	AGH143092	0.20 / 0.36 (0.008 / 0.014)	17	17	180	Straight
Grand Slam	AG141002	0.36 (0.014)	4	4	180	Straight
	AG141302				300	Straight
	AG141002J				180	J
ASAHI RG3	AHW10S302S	0.26 (0.010)	3	8	330	Straight
ASAHI EXTENSION 150cm	AG149000	0.36 (0.014)	—	—	150	—
ASAHI EXTENSION 165cm	AG149001	0.36 (0.014)	—	—	165	—



*ASAHI Conquest Pro 8-20 / ASAHI Confianza Pro 8-20 is not CE approved.

Not all products and shapes are available in all markets. Please check availability with your sales representative.
Coated with SLIP-COAT® coating.

© 2017 ASAHI INTECC CO., LTD.
"ASAHI", "ACT ONE", "ASAHI SION", "ASAHI Gaia", "Fielder", and "SUOH" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

Manufactured by:

Distributed by:



ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan
TEL:+81-561-48-5551 FAX:+81-561-48-5552

<http://www.asahi-intecc.com>