


SuperCross Microcatheters

Superb Crossability in Tortuous Anatomy

Superb Crossability in Challenging Cases

SuperCross Microcatheters provide support in tortuous anatomy and navigating bifurcated vessels.


Angled Tip Catheters

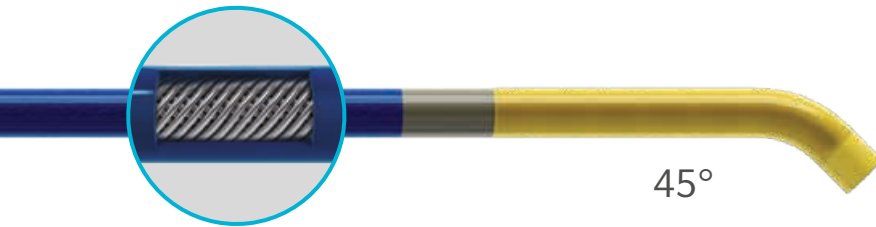

45°

Dual Coil Shaft Construction


Provides excellent torque response, flexibility, pushability and kink resistance

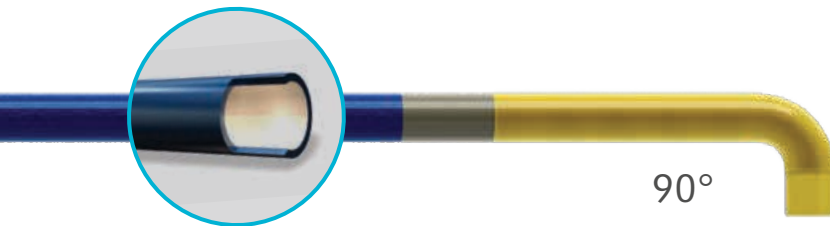

90°

PTFE Inner Layer

Allows for outstanding guidewire movement and delivery

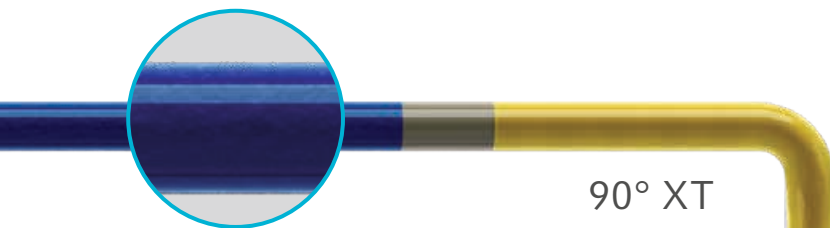

90° XT

XT=90° Extended Tip for
Secure Cannulation

Hydrophilic Coating on the Distal 80 cm

Enables exceptional trackability around tight bends and tortuous anatomy



120°

Embedded Platinum/Tungsten Coil

Provides enhanced visibility along the entire angled tip

Straight | Flexible Tip (FT) Catheters



Straight | FT

Straight Tip Version

Stiffer tip for more supportive delivery and increased pushability to access complex anatomy

Stainless Steel Braid

For flexibility, pushability and kink resistance

Flexible Tip (FT) Version

For increased trackability through tortuous anatomy


Ordering Information

SuperCross Microcatheters

The SuperCross Microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

MODEL	GUIDEWIRE COMPATIBILITY	DISTAL TIP	WORKING LENGTH	SHAFT CONSTRUCTION	DISTAL O.D.	PROXIMAL O.D.	DISTAL I.D.	PROXIMAL I.D.	HYDROPHILIC COATING
5302	0.014"	45°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5303	0.014"	45°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5304	0.014"	90°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5305	0.014"	90°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5306	0.014"	120°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5307	0.014"	120°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5308	0.014"	90° XT	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5309	0.014"	90° XT	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm

Packaged in quantities of 1 unit per box.

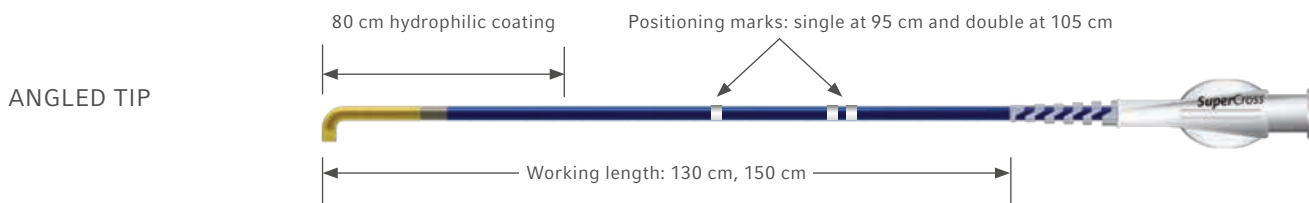


SuperCross Microcatheters—Straight | FT

The SuperCross Microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

MODEL	GUIDEWIRE COMPATIBILITY	DISTAL TIP	WORKING LENGTH	SHAFT CONSTRUCTION	DISTAL O.D.	PROXIMAL O.D.	DISTAL I.D.	PROXIMAL I.D.	HYDROPHILIC COATING
5300	0.014"	Straight	130 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm
5301	0.014"	Straight	150 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm
5340	0.014"	FT	130 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm
5341	0.014"	FT	150 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm

Packaged in quantities of 1 unit per box.

FT= Flexible Tip

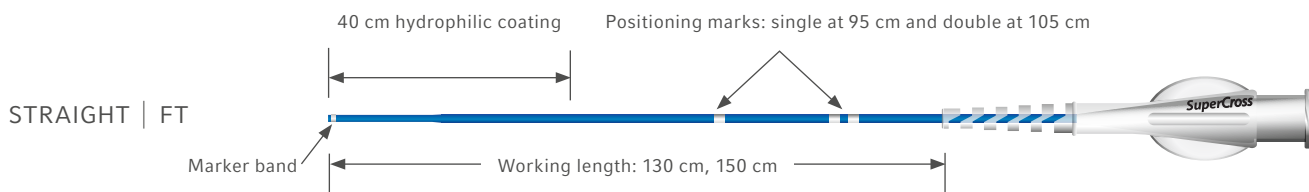

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose-driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit teleflex.com.

Teleflex is the home of Arrow, Deknatel, LMA, Pilling, QuikClot, Rüschi, UroLift, and Weck – trusted brands united by a common sense of purpose.

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Latin America: Phone +1 919 433 4999, la.cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

International: Phone +353 (0)9 06 46 08 00, orders.intl@teleflex.com, Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

Australia 1300 360 226

Austria +43 (0)1 402 47 72

Belgium +32 (0)2 333 24 60

Canada +1 (0)800 387 9699

China (Shanghai) +86 (0)21 6060 7000

China (Beijing) +86 (0)10 6418 5699

Czech Republic +420 (0)495 759 111

France +33 (0)5 62 18 79 40

Germany +49 (0)7151 691 3004

Greece +30 210 67 77 717

India +91 (0)80 4093 4790

Italy +39 0362 58911

Japan +81 (0)3 6632 3600

Korea +82 2 536 7550

Mexico +52 55 5002 3500

Netherlands +31 (0)88 00 215 00

New Zealand 0800 601 100

Poland +48 22 4624032

Portugal +351 225 491 051

Singapore (SEA non-direct sales countries) +65 6439 3000

Slovak Republic +421 (0)3377 254 28

South Africa +27 (0)11 807 4887

Spain +34 91 198 84 31

Switzerland +41 (0)31 818 40 90

United Kingdom +44 (0)1494 53 27 61

For more information, please visit teleflex.com

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

Teleflex, the Teleflex logo, Arrow, Deknatel, LMA, Pilling, QuikClot, Rüschi, SuperCross, UroLift and Weck are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 03/2022.

© 2022 Teleflex Incorporated. All rights reserved.

MCI-2018-0396-EN- REV 1 - 03 22 PDF